



CT Lumbar Spine With IV Contrast POST MYELOGRAM

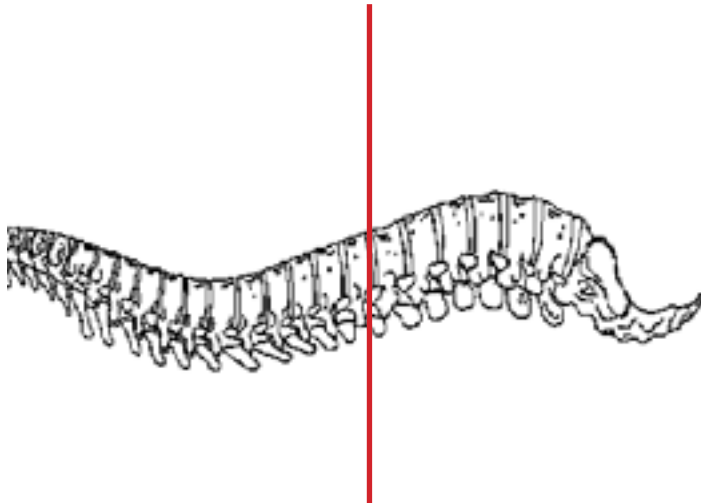

Pt Position	Supine
Contrast	Administered during myelogram
Injection Rate	N/A
Respiration	Quiet breathing
Volume Aquisition Specs	Appropriate to achieve images as specified in the following tables Small FOV T 11/12 to tip of coccyx

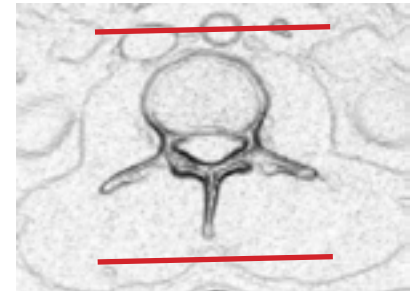

PLANE	ALGORITHM	**THICKNESS	DFOV	
Axial +C	Bone	2.5 mm	approx 180 mm	
Axial +C	Soft Tissue	2.5 mm	approx 180 mm	
Coronal +C	Bone	2.5 mm	approx 180 mm	reformat
Coronal +C	Soft Tissue	2.5 mm	approx 180 mm	reformat
Sagittal +C	Bone	2.5 mm	approx 180 mm	reformat
Sagittal +C	Soft Tissue	2.5 mm	approx 180 mm	reformat

** 2.5 mm slice thickness preferred / minimum 3 mm slice thickness