

## MRI BREAST IMPLANTS DYNAMIC


### Scan Range

FOV fit to patient

### Patient Position

Localizer 3 plane scout

Plane	Sequence	Fat Sat	Mode	SLT / SP	FOV	NS	GAP	NAQ
Axial	SILICONE SUP			4 mm	36x25	40	0.8	1
Axial	SALINE SUP			4 mm	36x25	79	0.8	1
Axial	PRE DYNAMIC			2.4 mm	34X26	79		1
Axial	DYNAMIC			2.4 mm	34x26	79		1
Sagittal	HIRES			3 mm	25	83		1

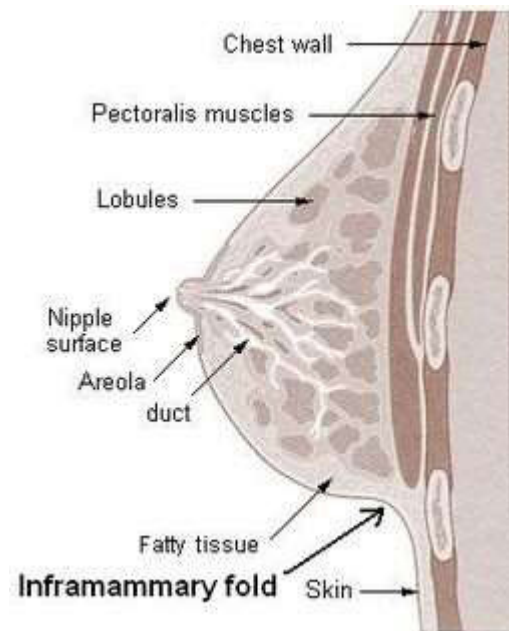

### SPECIAL INSTRUCTIONS

- Menstrual cycle 5-15 days for routine patients or follow-ups
- Positioning is KEY to breast exam dictation. Ensure there is no inframammary fold at the bottom of the breast. No extreme compression, sharp edges, or folds in breast tissue.
- **NO BELLY SKIN BULGE – PULL BREAST STRAIGHT FOR 90 DEGREE ANGLE**