


# MRI CHEST MEDIASTINAL MASS WITHOUT & WITH CONTRAST

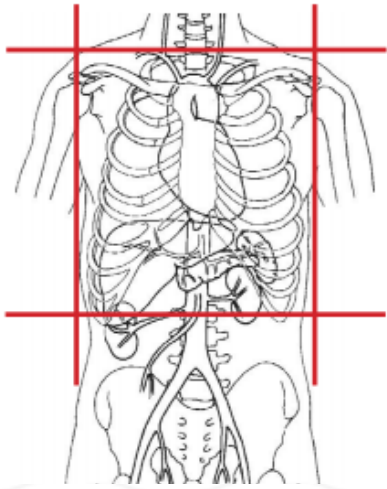

Indications	Mediastinal mass
Patient Position	Supine, arms out of FOV
Contrast	Without and with
Respiration	Breath- hold should be encouraged
Scan Parameters	Cover mediastinum and mass
Special notes	If unsure of scan parameters, check with Radiologist.

PLANE	Sequence	Fat Sat	Mode	SLT/SP	FOV	NOTES
Sagittal	T2	Y			fit to patient	
Sagittal	T1				fit to patient	
Axial	T1				fit to patient	
Axial	T2				fit to patient	
Axial	T2	Y			fit to patient	
Axial	3D T1		FSPGR VIBE LAVA			Ultra fast GRE/ your scanner equivalent
Sagittal	3D T1		FSPGR VIBE LAVA			Ultra fast GRE/ your scanner equivalent
Sagittal + C	3D T1		FSPGR VIBE LAVA			Ultra fast GRE/ your scanner equivalent
Axial + C	3D T1		FSPGR VIBE LAVA			Ultra fast GRE/ your scanner equivalent

Any deviation from protocol MUST be Radiologist approved. Rad Director will be notified if this occurs without prior approval.