


MRI SACRUM WITHOUT & WITH CONTRAST

Indications	Sacral fracture, sacral mass, SI joint infection, NOT for female reproductive or ovarian complaints.
Patient Position	Supine- feet internally rotated and patient made comfortable with sponges. Do not rotate if suspicion of femoral neck / hip fx.
FOV	Fits to patient/ body part
Localizer	3 plane scout localize
Notes	Include both sides of sacrum for the Sagittal sequences

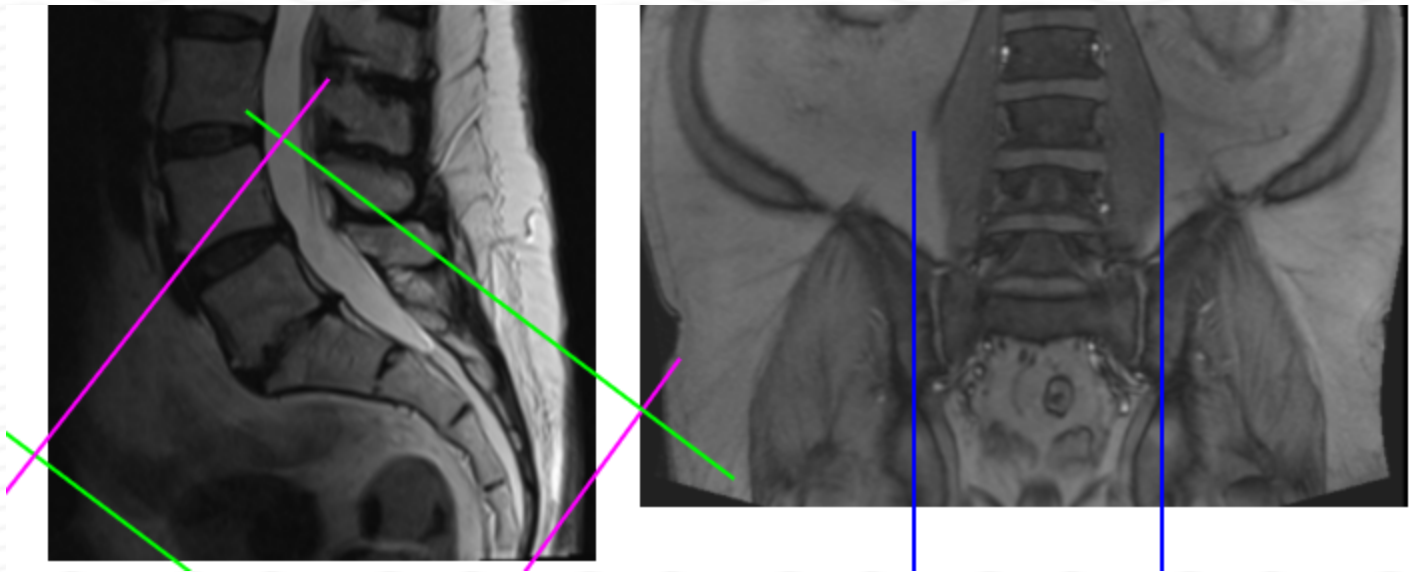

PLANE	SEQUENCE	FAT SAT	MODE	SLT/SP	FOV	NOTES
Sagittal	T1 DIXON TSE			4 / 1	Approx 220 mm	
Axial oblique	T1 Dixon TSE			4 / 1	Approx 180 mm	
Axial oblique	T1	Y		4 / 1	Approx 180 mm	
Sagittal	STIR TSE			4 / 1	Approx 220 mm	
Coronal oblique	STIR TSE			4 / 1	Approx 180 mm	
INJECT CONTRAST				4 / 1		
Sagittal + C	T1 Dixon TSE	Y		4 / 1	Approx 220 mm	
Axial oblique + C	T1 Dixon TSE			4 / 1	Approx 180 mm	