


## **CT BONY PELVIS without IV CONTRAST**

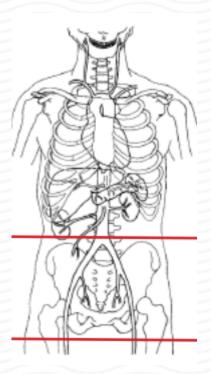

Indications	possible fracture
Contraindications	Pregnancy is a relative contraindication
Patient Position	Supine
Contrast	NONE
Respiration	N/A
Scan Parameters	Iliac crest through lesser trochanter

PLANE	ALGORITHM	THICKNESS	DFOV	
Axial	Bone	3.0 mm	pt largest + 4 cm	
Axial	Soft tissue	3.0 mm		
Coronal & Sagittal	Bone	3.0 mm		Reformats

## NOTES:

• LMP on pts of child-bearing age

· Shield children when possible

Any deviation from protocol MUST be Radiologist approved. Rad Director will be notified if this occurs without prior approval.