
Reviewed: Jan 2024 Updated: Sept 2020



CT BONY PELVIS without IV CONTRAST

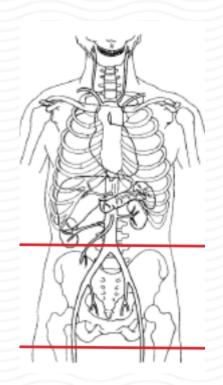

Indications	possible fracture
Contraindications	Pregnancy is a relative contraindication
Patient Position	Supine
Contrast	NONE
Respiration	N/A
Scan Parameters	Iliac crest through lesser trochanter

PLANE	ALGORITHM	THICKNESS	DFOV	
Axial	Bone	3.0 mm	pt largest + 4 cm	
Axial	Soft tissue	3.0 mm		
Coronal & Sagittal	Bone	3.0 mm		Reformats

NOTES:

- LMP on pts of child-bearing age
- Shield children when possible