
Reviewed: Jan 2024 Revised: Mar 2022 **New: Jan 2022**



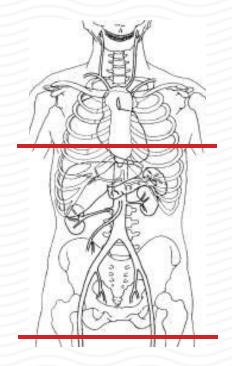

CTA RENAL ARTERY

Indications	evaluation of renovascular hypertension, presurgical planning , vascular interventional planning		
Patient Position	Supine		
Contrast	Facility specific based on weight/ renal function		
Injection rate	20 g or larger in large vein (AC). 4 ml per sec, 120 mL based on weight / renal function. Prefer 18 g		
Respiration	Breath- hold		
Scan Parameters	Diaphragm to bottom of femoral heads.		

PLANE	ALGORITHM	THICKNESS	DFOV	
Axial + C	Abdomen	3.0 mm		EARLY ARTERIAL BOLUS TRACKING
Coronal and Sagittal + C	Abdomen	2 mm x 2 mm		Reformats of contrast run above
Axial, Coronal + C	MIP	7.0 mm x 2.00 mm		MIPS of early arterial run.
Axial + C	Abdomen	3.0 mm x 3.00 mm		90 SEC DELAY