


MRI BREAST DYNAMIC

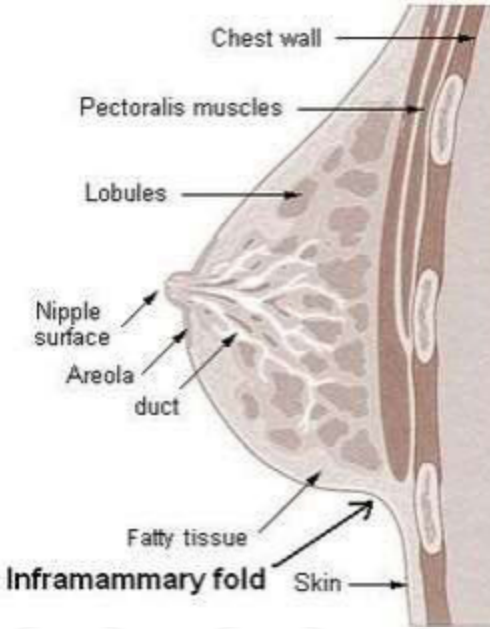

| | |
|------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Scan Parameters | Fit to patient |
| Special notes | Menstrual cycle 5-15 days for routine patients or follow-ups |
| | Positioning is KEY to breast exam dictation. Ensure there is no inframammary fold at the bottom of the breast. No extreme compression, sharp edges, or folds in the breast tissue. |
| | No belly skin bulge- pull breast straight for 90 degree angle. |

| PLANE | Sequence | Fat Sat | Mode | SLT/SP | FOV | NS | GAP | NAQ |
|----------|----------|---------|------|--------|---------|----|-----|-----|
| Axial | T1 | No | | 4 mm | 37 x 24 | 40 | 0.8 | 1 |
| Axial | T2 | Yes | | 4 mm | 37 x 24 | 40 | 0.8 | 1 |
| Axial | Pre DYN | | | 2.4 mm | 34 x 26 | 79 | | 1 |
| Axial | DYN | | | 2.4 mm | 34 x 26 | 79 | | 1 |
| Sagittal | HIRES | | | 3 mm | 25 | 83 | | 1 |