


MRI NECK WITHOUT

Use this protocol for the following:	THYROID/CERVICAL TRACHEA AND CERVICAL ESOPHAGUS
Scan Range	Variable according to area of interest/ body part.
FOV	Fits to patient/ body part
Notes	Use skin marker over area of tenderness if no palpable mass and above / below area of mass

PLANE	SEQUENCE	FAT SAT	MODE	SLT/SP	FOV	NOTES
Axial	T1	No		3.0 / 0	240 mm	
Axial	T2	No		3.0 / 0	240 mm	
Coronal angled	T2	Yes		3.0 / 0	240 mm	
Axial	DWI	No		3.0 / 0	240 mm	
Sagittal	T2	No		3.0 / 0	240 mm	
Axial	T2	Yes		3.0 / 0	240 mm	

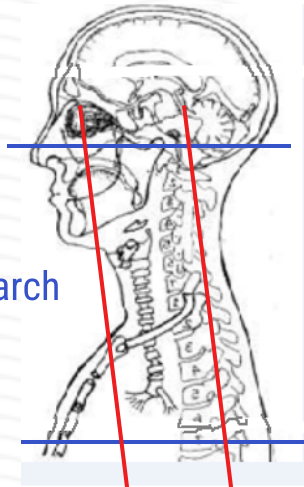

Hard palate to below aortic arch

