

Reviewed: Jan 2025 Revised: Sept 2020



CT SINUS / FACIAL BONES WITHOUT IV CONTRAST



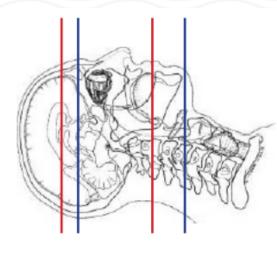

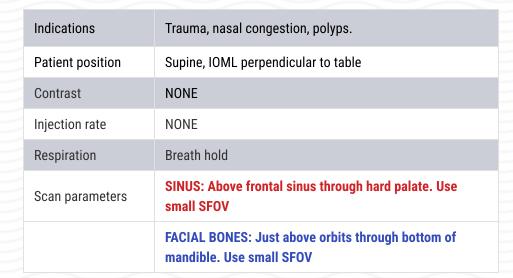



PLANE	ALGORITHM	THICKNESS	DFOV	
Axial	Soft tissue	2.0 mm	170- 190 mm	
Axial	Bone	2.0 mm		
Coronal	Bone	2.0 mm		Reformats: Frontal sinus through mastoids
Sagittal	Bone	2.0 mm		Reformats: Cover all bony structures

NOTES:

- · LMP on pts of child-bearing age
- Shield children when possible
- · Follow GFR and metformin protocol for contrast studies

Any deviation from protocol MUST be Radiologist approved. Rad Director will be notified if this occurs without prior approval.