


MRI CHEST - MEDIASTINAL MASS WITHOUT & WITH CONTRAST

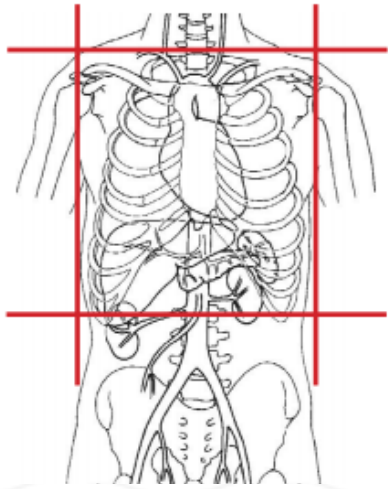

Indications	Mediastinal mass
Patient Position	Supine, arms out of FOV
Contrast	Without and with
Respiration	Breath- hold should be encouraged
Scan Parameters	Cover mediastinum and mass
Special notes	If unsure of scan parameters, check with Radiologist.

PLANE	Sequence	Fat Sat	Mode	FOV	NOTES
Sagittal	T2	Y		fit to patient	
Sagittal	T1			fit to patient	
Axial	T1			fit to patient	
Axial	T2			fit to patient	
Axial	T2	Y		fit to patient	
Axial	3D T1		FSPGR VIBE LAVA		Ultra fast GRE/ your scanner equivalent
Sagittal	3D T1		FSPGR VIBE LAVA		Ultra fast GRE/ your scanner equivalent
Sagittal + C	3D T1		FSPGR VIBE LAVA		Ultra fast GRE/ your scanner equivalent
Axial + C	3D T1		FSPGR VIBE LAVA		Ultra fast GRE/ your scanner equivalent